 <p>B.S. Aquino Drive, Bacolod City, Negros Occidental, 6100</p> <p>DR. PABLO O. TORRE MEMORIAL HOSPITAL</p>	Document Code:	DPOTMH-E-66-P02-S03
	Effective Date:	06-30-2022
	Document Type:	Standard Operating Procedure
	Page Number:	1 of 3
	Department/Section:	Physical Medicine & Fitness Center
	Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

### **PURPOSE:**

To describe in detail the steps in the operation of the gym equipment to ensure accuracy and safety while using the machine.

### **SCOPE:**


Applies to all Physical Medicine & Fitness Center Staff of Dr. Pablo O. Torre Memorial Hospital (DPOTMH)

### **RESPONSIBLE PERSON:**


Physical Therapist, Fitness Instructor

### **PROCEDURE:**


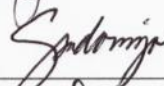
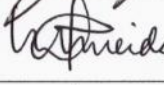





1. Physiatrist clears the client as candidate for a fitness program.
2. PT/FI will assess the client's fitness level and will formulate the exercise program.
3. FI/PT checks the client's BP and HR if necessary prior to the start of exercise.
4. To begin a work out, grasp the handles and carefully step on the pedals.
5. Keep the shoes secured by the bike pedal strap during a work out.
6. To begin a work out, press the START key or begin pedaling to activate the console for programming a new work out. The activated console lights up and displays SELECT WORKOUT OR PRESS QUICK START. If no key is pressed for 15 seconds, the console shuts down again.
7. Press the ENTER key after entering each value, such as WEIGHT or WORK OUT TIME, when prompted by the console display messages.
8. Use the numeric keypad to enter values such as LENGTH OF WORK OUT, AGE, TARGET HEART RATE, and INTENSITY LEVEL. While a work out is in progress, numeric keys can be used to set new Intensity Level or change the target heart rate.
9. Use the ARROW keys for the similar function.
10. Press QUICK START to begin a manual work out immediately.
11. Press CLEAR key to clear incorrect data. Pressing CLEAR the second time causes the MESSAGE CENTER to prompt the user for a new work out set – up.

 <p>B.S. Aquino Drive, Bacolod City, Negros Occidental, 6100</p> <p>DR. PABLO O. TORRE MEMORIAL HOSPITAL</p>	Document Code:	DPOTMH-E-66-P02-S03
	Effective Date:	06-30-2022
	Document Type:	Standard Operating Procedure
	Page Number:	2 of 3
	Department/Section:	Physical Medicine & Fitness Center
	Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

12. Press the following work out set – up as desired: FAT BURN, CARDIO, RANDOM, MANUAL, ZONE TRAINING, PERSONAL TRAINER, and HILL PLUS.
13. Work out program ends automatically in a certain COOL DOWN MODE, which lowers the intensity level. Press the COOL DOWN key at any desired point during a work out to go immediately into COOL DOWN MODE.
14. The MESSAGE CENTER window displays step – by – step instructions for setting up a work out. It displays statistics about the progress of the work out: HR, Distance, Speed, and Calories.
15. During a work out, grasp two of the tour sensors that are not in the handle bars to HR monitoring. Maintain a vertical trunk alignment as much as possible.
16. Begin the pedaling when time starts counting the TIME WINDOW.
17. The exercise ends when COOL DOWN time elapsed.
18. WORK OUT SUMMARY appears at the MESSAGE CENTER at the end of the exercise.
19. If, at any time while exercising, the user experiences faintness, dizziness, pain or shortness of breath, he or she must stop immediately.

 <p>B.S. Aquino Drive, Bacolod City, Negros Occidental, 6100</p> <p>DR. PABLO O. TORRE MEMORIAL HOSPITAL</p>	Document Code:	DPOTMH-E-66-P02-S03
	Effective Date:	06-30-2022
	Document Type:	Standard Operating Procedure
	Page Number:	3 of 3
	Department/Section:	Physical Medicine & Fitness Center
	Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

### APPROVAL:

	Name/Title	Signature	Date
Prepared by:	<b>MA. CLAUDETTE L. SOLANCHO, PTRP</b> PM&FC Staff		7/11/22
	<b>SUZETTE M. DOMINGO, PTRP</b> PM&FC Staff		7-11-22
Verified:	<b>MA. TERESA O. ANIERDES, PTRP, MBA-HA</b> PM&FC Supervisor		7/11/22
	<b>JOHN ANDREW B. SANCHEZ, MD, FPARM</b> Physiatrist		7/11/22
Reviewed by:	<b>DENNIS C. ESCALONA, RN, MN, FPSQua</b> Quality Assurance Supervisor		7/11/2022
Recommending Approval:	<b>ROSARIO D. ABARING, RN, MAN, PhD</b> Ancillary Services Division Officer		07.11.2022
	<b>FREDERIC IVAN L. TING, MD</b> OIC-Total Quality Division		7/18/22
Final approved by:	<b>GENESIS GOLDI D. GOLINGAN</b> President & CEO		9/29/22






DR. PABLO O. TORRE  
MEMORIAL HOSPITAL

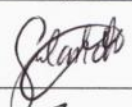
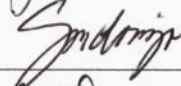
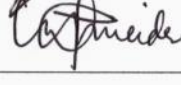





B.S. Aquino Drive,  
Bacolod City,  
Negros Occidental,  
6100

Document Code:	DPOTMH-E-66-P02-WI03
Effective Date:	06-30-2022
Document Type:	Work Instruction
Page Number:	1 of 2
Department/Section:	Physical Medicine & Fitness Center
Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

KEY TASKS	PERSON RESPONSIBLE
1. Clears the client as candidate for a fitness program.	Physiatrist
2. Assesses the client's fitness level and formulates the exercise program.	
3. Grasps the handles to begin a work out and carefully steps on the pedals.	Physical Therapist/Fitness Instructor
4. Keeps the shoes secured by the bike pedal strap during a work out.	
5. Begins the work out by pressing the START key or begin pedaling to activate the console for programming a new work out.	
6. Enters values such as LENGTH OF WORK OUT, AGE, TARGET HEART RATE, and INTENSITY LEVEL.	
7. Presses the following work out set-up as desired: FAT BURN, CARDIO, RANDOM, MANUAL, ZONE TRAINING, PERSONAL TRAINER, and HILL PLUS.	Client
8. Presses the COOL DOWN key at any desired point during a work out to go immediately into COOL DOWN MODE.	
9. Advises client to stop immediately if he/she experiences faintness, dizziness, pain or shortness of breath at any time while exercising.	

 <p>DR. PABLO O. TORRE MEMORIAL HOSPITAL</p> <p>B.S. Aquino Drive, Bacolod City, Negros Occidental, 6100</p>	Document Code:	DPOTMH-E-66-P02-WI03
	Effective Date:	06-30-2022
	Document Type:	Work Instruction
	Page Number:	2 of 2
	Department/Section:	Physical Medicine & Fitness Center
	Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

**APPROVAL:**

	Name/Title	Signature	Date
Prepared by:	<b>MA. CLAUDETTE L. SOLANCHO, PTRP</b> PM&FC Staff		7/11/22
	<b>SUZETTE M. DOMINGO, PTRP</b> PM&FC Staff		7-11-22
Verified:	<b>MA. TERESA O. ANIERDES, PTRP, MBA-HA</b> PM&FC Supervisor		7/11/22
	<b>JOHN ANDREW B. SANCHEZ, MD, FPARM</b> Physiatrist		7/11/22
Reviewed by:	<b>DENNIS C. ESCALONA, RN, MN, FPSQua</b> Quality Assurance Supervisor		7/11/2022
Recommending Approval:	<b>ROSARIO D. ABARING, RN, MAN, PhD</b> Ancillary Services Division Officer		07.11.2022
	<b>FREDERIC IVAN L. TING, MD</b> OIC-Total Quality Division		7/18/22
Final approved by:	<b>GENESIS GOLDI D. GOLINGAN</b> President & CEO		7/29/22

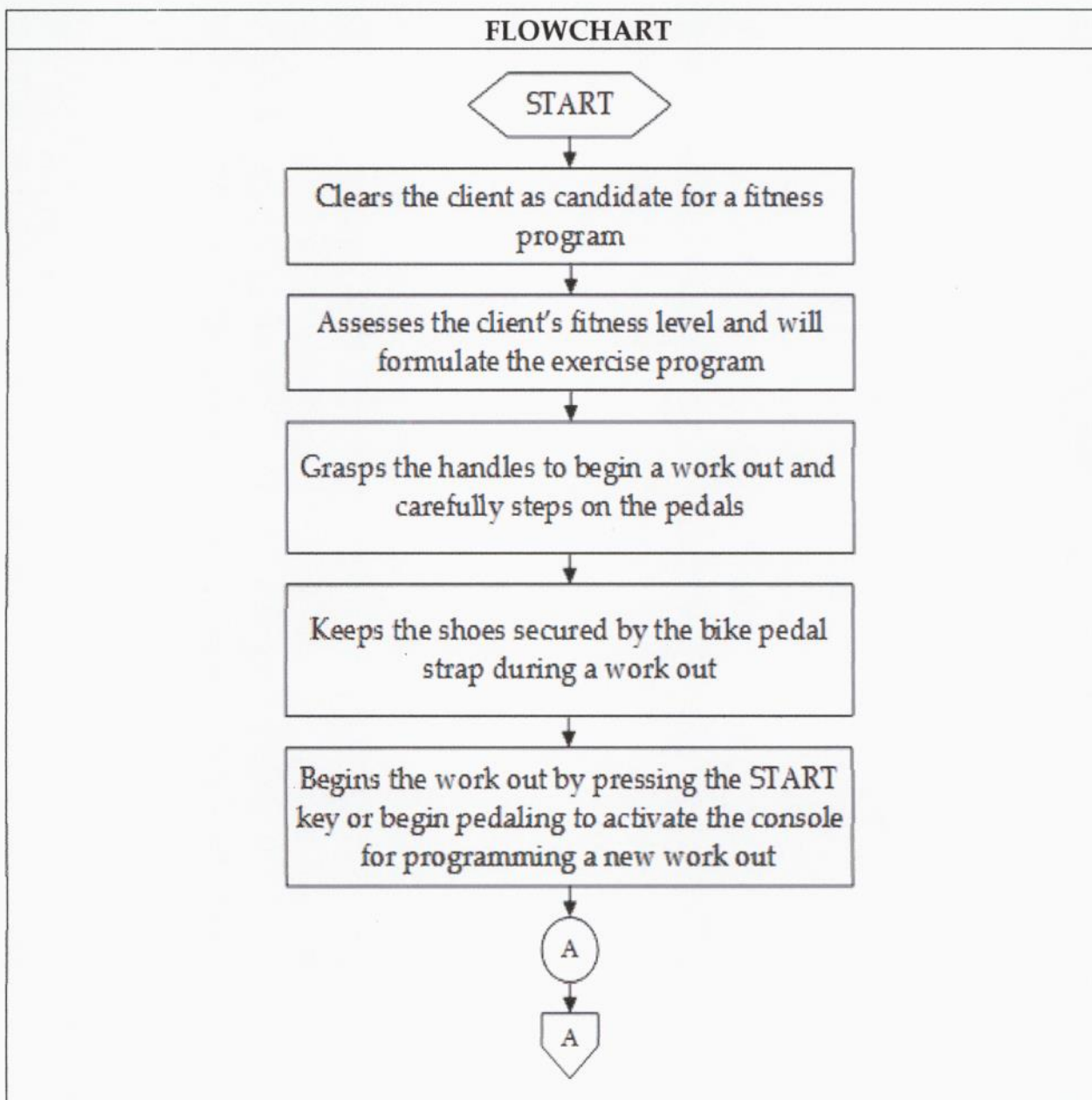


DR. PABLO O. TORRE  
MEMORIAL HOSPITAL

B.S. Aquino Drive,  
Bacolod City,  
Negros Occidental,  
6100

Document Code:	DPOTMH-E-66-P02-FC03
Effective Date:	06-30-2022
Document Type:	Flowchart
Page Number:	1 of 3
Department/Section:	Physical Medicine & Fitness Center
Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

## FLOWCHART





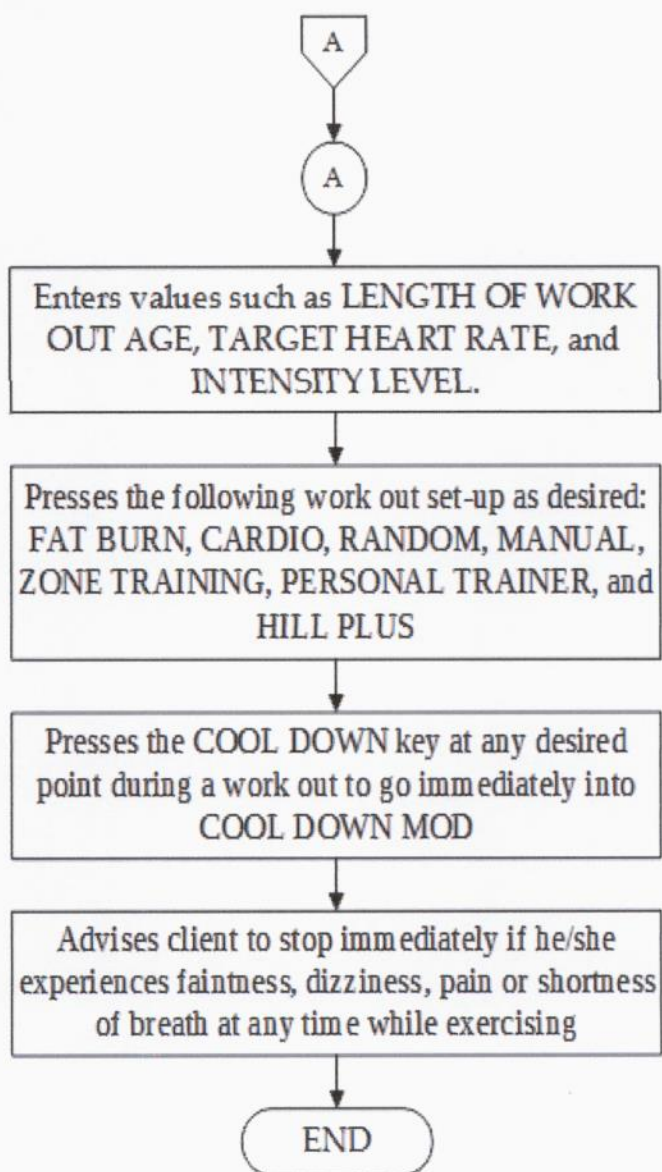



DR. PABLO O. TORRE  
MEMORIAL HOSPITAL

B.S. Aquino Drive,  
Bacolod City,  
Negros Occidental,  
6100

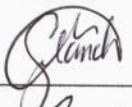

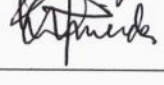


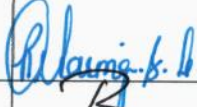


Document Code:	DPOTMH-E-66-P02-FC03
Effective Date:	06-30-2022
Document Type:	Flowchart
Page Number:	2 of 3
Department/Section:	Physical Medicine & Fitness Center
Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

### FLOWCHART



 <p><b>DR. PABLO O. TORRE MEMORIAL HOSPITAL</b></p> <p>B.S. Aquino Drive, Bacolod City, Negros Occidental, 6100</p>	Document Code:	DPOTMH-E-66-P02-FC03
	Effective Date:	06-30-2022
	Document Type:	Flowchart
	Page Number:	3 of 3
	Department/Section:	Physical Medicine & Fitness Center
	Document Title:	<b>GYM EQUIPMENT PROTOCOL: (95 xi TOTAL BODY TRAINER)</b>

**APPROVAL:**

	Name/Title	Signature	Date
Prepared by:	<b>MA. CLAUDETTE L. SOLANCHO, PTRP</b> PM&FC Staff		7/11/22
	<b>SUZETTE M. DOMINGO, PTRP</b> PM&FC Staff		7-11-22
Verified:	<b>MA. TERESA O. ANIERDES, PTRP, MBA-HA</b> PM&FC Supervisor		7/11/22
	<b>JOHN ANDREW B. SANCHEZ, MD, FPARM</b> Physiatrist		7/11/22
Reviewed by:	<b>DENNIS C. ESCALONA, RN, MN, FPSQua</b> Quality Assurance Supervisor		7/11/2022
Recommending Approval:	<b>ROSARIO D. ABARING, RN, MAN, PhD</b> Ancillary Services Division Officer		07.11.2022
	<b>FREDERIC IVAN L. TING, MD</b> OIC-Total Quality Division		7/18/22
Final approved by:	<b>GENESIS GOLDI D. GOLINGAN</b> President & CEO		9/29/22